



ELSEVIER

Contents lists available at ScienceDirect

MethodsX


journal homepage: www.elsevier.com/locate/mex

Step-by-step quantitative analysis of focal adhesions

Utku Horzum, Berrin Ozdil, Devrim Pesen-Okvur*

Izmir Institute of Technology, Department of Molecular Biology and Genetics, 35430 Izmir, Turkey

GRAPHICAL ABSTRACT

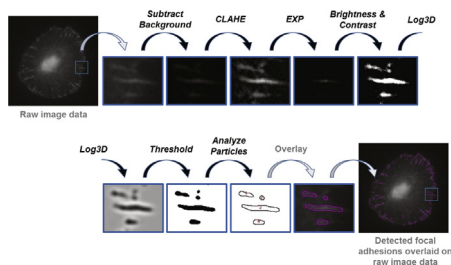

ABSTRACT

Focal adhesions (FAs) are specialized adhesive structures which serve as cellular communication units between cells and the surrounding extracellular matrix. FAs are involved in signal transduction and actin cytoskeleton organization. FAs mediate cell adhesion, which is a critical phenomenon in cancer research. Since cells can form many and micrometer scale FAs, their quantitative analysis demands well-optimized image analysis approaches [1–3]. Here, we have optimized the analysis of FAs of MDA-MB-231 breast cancer cells. The optimization is based on proper processing of immunofluorescence images of vinculin, which is one of the markers of FAs. All image processing steps are carried out using the ImageJ software, which is freely available and in the public domain. The advantages of our method are:

- The analysis steps are simplified by combining different plugins of the ImageJ program.
- FAs are better detected with minimal false negatives due to optimized processing of fluorescent images.
- This approach can be applied to quantify a variety of fluorescent images comprising focal and/or localized signals within a high background such as FAs, one of the many complex signaling structures in a cell.

© 2014 The Authors. Published by Elsevier B.V. This is an open access article under the CC BY license (<http://creativecommons.org/licenses/by/3.0/>).

* Corresponding author.

E-mail addresses: devrimpesen@iyte.edu.tr, devrimpesenokvur@gmail.com (D. Pesen-Okvur).

ARTICLE INFO

Method name: Analysis of fluorescent images of focal adhesions

Keywords: Focal adhesion, Particle analysis, Image processing, ImageJ

Article history: Received 21 March 2014; Accepted 27 June 2014

Materials

- Fluorescence images
- ImageJ software

Critical points

- Sample should be properly immunostained so that good quality fluorescence images can be captured.
- The images should be of high quality; unfocused images and inadequate resolution hinder the quality of the analysis.
- Fluorescence images should ideally be captured at 100× magnification.
- If color images are captured, single channel images should first be created with the ImageJ `SPLIT CHANNEL` command.
- Use raw image data which are greyscale images without any special preprocessing.

Method details

Here, we have analyzed focal adhesions of MDA-MB-231 cells. However, our method can be applied to not only any cell type such as fibroblasts, endothelial cells or leucocytes but also to any structure in the cell with a localized signal such as focal adhesions, endosomes or invadopodia. Focal adhesions were immunostained for vinculin using a vinculin specific primary antibody (Sigma Cat no: V9131), followed by a fluorophore conjugated secondary antibody. Other markers of focal adhesions such as paxillin or focal adhesion kinase could also be easily used for analysis. 8-bit images were captured using an Olympus epi-fluorescence microscope with an Olympus infinity corrected 100× oil UPlanSApo objective with a numerical aperture of 1.4 and an Olympus SLR E-330 camera such that each image was 3136 pixels by 2352 pixels and 37 pixels corresponded to 2 μm. All steps of image processing were carried out using ImageJ which is freely available and in the public domain. ImageJ provides a wide range of processing and analysis approaches possible via not only built-in functions but also numerous plugins. ImageJ has been in use for some time by researchers. However, choosing the right combination and order of processing steps is crucial for obtaining meaningful quantitative results.

The image processing optimized here aims to identify FAs immunostained for vinculin. The raw fluorescent images are processed in steps as follows:

- Step 1: Apply `SUBTRACT BACKGROUND`. Selection of sliding paraboloid option changes the rolling ball to a parabolic with the same ball radius in pixels but with sharper curvature. A parabola slides in different directions over the image and calculates and subtracts local background from the original image. We chose the `SLIDING PARABOLOID` option with the `ROLLING BALL` radius set to 50 pixels.
- Step 2: Enhance the local contrast of the image by running `CLAHE` (Contrast Limited Adaptive Histogram Equalization). The plugin `CLAHE` offers three critical parameters which are block size, histogram bins and maximum slope to clarify desired objects. We used the following values: block size=19, histogram bins=256, maximum slope=6, no mask and fast [4].
- Step 3: Apply mathematical exponential (`EXP`) to further minimize the background.
- Step 4: Adjust `BRIGHTNESS & CONTRAST` automatically. The `BRIGHTNESS & CONTRAST` tool updates the lookup table of the image based on an analysis of the histogram of the image.

- Step 5: Run `Log3D` (Laplacian of Gaussian or Mexican Hat) filter. The `Log3D` plugin filters the image based on user predefined parameters which are standard deviations in X, Y and Z directions. For 2D images Z direction is not relevant and slice per slice processing is used. We defined the size of `Log3D` filter as sigma X=5 and sigma Y=5 [5].
- Step 6: Run the `THRESHOLD` command. The `THRESHOLD` command converts a grayscale image to a binary image with two pixel values, 255 (white) and 0 (black). For the threshold command, there are options such as Huang, triangle, mean or default. In addition, the user can adjust minimum and maximum threshold values or set the threshold levels according to an analysis of the histogram of the current image automatically. Here we used the default method and threshold was adjusted automatically.
- Step 7: Execute `ANALYZE PARTICLES` command. The `ANALYZE PARTICLES` command scans the thresholded (binary) image and finds the edges of objects or particles. Essentially, particle analyzer counts and measures particles according to user pre-defined parameters which are size and circularity. We set these parameters as follows: size=50-infinity and circularity=0.00–0.99.

As a result of the analyze particles process, outlines of the particles found can be saved and overlaid on the original image to visually confirm faithful detection of FAs.

Our analysis on FAs of MDA-MB-231 cells cultured on fibronectin, an extracellular matrix protein that promotes cell adhesion and on K-casein, a blocking protein that prevents cell adhesion, determined the following [6]:

	Number of FAs per cell	Area of FA	Total area of FAs per cell
Fibronectin	65 ± 14	0.56 ± 0.02 μm ²	36.09 ± 8.45 μm ²
K-casein	5 ± 1	0.38 ± 0.04 μm ²	2.04 ± 0.63 μm ²

The data is presented as average ± standard error.
Number of FAs per cell=Total number of FAs identified in one cell.
Area of FA=Area of a single FA.
Total area of FAs per cell=The area of all FAs in one cell are summed up.
Number of cells analyzed=6–9 number of FAs analyzed=0–348.

As expected, cells form more and larger FAs on a fibronectin coated surface compared with a K-casein coated surface ($p < 0.05$).

Additional information

Background information

Focal adhesions (FAs) are multifunctional structures comprising over a hundred proteins. FAs provide not only physical contact between a cell and its extracellular environment but also mediate inside-out and outside-in signaling [7]. Therefore, quantitative analysis of FAs is crucial for understanding important cellular phenomena such as cell morphology, motility, tissue organization, growth and differentiation [8,9]. Since cells can form many and micrometer scale FAs, their quantitative analysis demands well-optimized image analysis approaches [1–3]. The dimensions of focal adhesions range from 0.25 to 10 μm depending on the maturation stage of the FA as well as the cell type. Mature focal adhesions are also associated with the cytoskeleton [2,10]. There can be more than a hundred proteins associated with a focal adhesion and the identification of focal adhesion is possible by labeling protein of interest within this complex. Labeling of a mechanical stabilizer protein, vinculin, is the best suitable method technically and commercially [10].

Useful links

ImageJ software: <http://rsbweb.nih.gov/ij/download.html>
Log3D plugin: <http://bigwww.epfl.ch/sage/soft/LoG3D/>
CLAE plugin: <http://rsbweb.nih.gov/ij/plugins/clahe/index.html>
Similar image processing package: <http://fiji.sc/Fiji>

Acknowledgements

This work was supported by TUBITAK Grant 111T026. *MethodsX* thanks the reviewers of this article for taking the time to provide valuable feedback.

References

- [1] M.E. Berginski, E.A. Vitriol, K.M. Hahn, S.M. Gomez, High-resolution quantification of focal adhesion spatiotemporal dynamics in living cells, *PLoS One* 6 (7) (2011) e22025, PubMed PMID: 21779367. Pubmed Central PMCID: Pmc3136503 (Epub 2011/07/23).
- [2] E. Lamers, R. van Horssen, J. te Riet, F.C. van Delft, R. Luttge, X.F. Walboomers, et al., The influence of nanoscale topographical cues on initial osteoblast morphology and migration, *Eur. Cells Mater.* 20 (2010) 329–343, PubMed PMID: 21061239 (Epub 2010/11/10).
- [3] E. Zamir, B.Z. Katz, S. Aota, K.M. Yamada, B. Geiger, Z. Kam, Molecular diversity of cell-matrix adhesions, *J. Cell Sci.* 112 (June (Pt 11)) (1999) 1655–1669, PubMed PMID: 10318759 (Epub 1999/05/13).
- [4] S. Saalfeld, CLAHE (Contrast Limited Adaptive Histogram Equalization), 2009 Available from: <http://rsbweb.nih.gov/ij/plugins/clahe/index.html> (updated 2009/11/17).
- [5] D. Sage, F.R. Neumann, F. Hediger, S.M. Gasser, M. Unser, Automatic tracking of individual fluorescence particles: application to the study of chromosome dynamics, *IEEE Trans. Image Process.* 14 (September (9)) (2005) 1372–1383, PubMed PMID: 16190472 (Epub 2005/09/30).
- [6] U. Horzum, B. Ozdil, D. Pesen-Okvur, Micrometer scale spacings between fibronectin nanodots regulate cell morphology and focal adhesions, *Mater. Res. Express* 1 (2014) 025402.
- [7] B. Geiger, K.M. Yamada, Molecular architecture and function of matrix adhesions, *Cold Spring Harb. Perspect. Biol.* 3 (5) (2011), PubMed PMID: 21441590. Pubmed Central PMCID: Pmc3101841 (Epub 2011/03/29).
- [8] A.L. Berrier, K.M. Yamada, Cell-matrix adhesion, *J. Cell. Physiol.* 213 (December (3)) (2007) 565–573, PubMed PMID: 17680633 (Epub 2007/08/08).
- [9] K. Burridge, M. Chrzanowska-Wodnicka, Focal adhesions, contractility, and signaling, *Annu. Rev. Cell Dev. Biol.* 12 (1996) 463–518, PubMed PMID: 8970735 (Epub 1996/01/01).
- [10] G.R. Owen, D.O. Meredith, I. ap Gwynn, R.G. Richards, Focal adhesion quantification – a new assay of material biocompatibility? Review, *Eur. Cells Mater.* 9 (2005) 85–96, discussion 85–96. PubMed PMID: 15977138 (Epub 2005/06/25).